

The treatment centre range to

suit all tastes

for better daily dentistry



Dental Chair/Package Range

8	EURUS S8 (the ultimate chair solution!)
20	CLESTA EIII SWING (ambidextrous and economical)
28	TBCOMPASS (ambidextrous with more features)
32	CLESTA elli
38	CLESTA II
44	EURUS Swivel chair (kneebreak)
47	CLAIR (kneebreak)

2	100 years of Belmont
48	Upholstery
50	Dental Light
54	Dental Xray
60	Stool
64	Dealers - Contact



Rise into beauty. Walk into life 100 years of Takara Belmont

TAKARA BELMONT have been contributing ourselves for people's wish to be healthy and beautiful through a century of our history.

Belmont Values

The principle of our craftsmanship is to realize 3 Cs for dentists

Cure

To treat patients with reliable technique and state-of-the-art technology

- Belmont's original electro-hydraulic technology achieves dental equipments with long-term & high durability.
- Ergonomically designed dental equipments make dentist's long working hours easy with minimal strain.
- Our dental equipments provide solutions against infection to stay clean and safe during treatment.





Care

To support preventative care and aesthetic dentistry

- It came to be widely known that the quality of life to be able to raise by keeping mouth healthy.
- For that purpose, a care and a maintenance are also important not a cure of a cavity.
- Our goal is to lessen the pressures felt during treatment, so that our dental equipment eliminates the noise, vibration and shock that lie at the root of feelings of insecurity among patients.

Communication

To enhance patient's consciousness and mutual understanding

- Patients can discuss needs and a thought with a dentist further in the relaxed space.
- With the aim of alleviating anxiety for patients visiting clinics and enhancing their relationships of trust with dentists, we have pursued dialog and comfort from the "patient's perspective" by reconfiguring the unit with the intent of "disguising" dental treatment.
- This "disguise" covertly stimulates the patient's senses in a good way.




Made by Hand and Mind, in Japan

Almost one hundred years on, we manufacture every piece of every chair, unit, light and x-ray unit we deliver to market.

The provenance of our product is thus unquestionable, and the high quality of our products are maintained by the familiarity of eye, hand and mind that comes with our loyal workers.





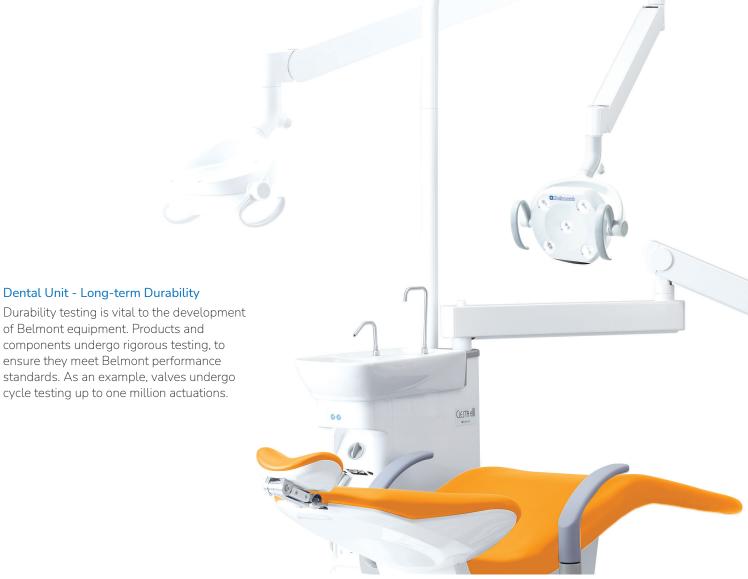



Dental Chair -Quality Starts Within

Each Belmont chair is the result of more than 50 years of product research and development, meticulous engineering expertise and innovative production methods. Belmont hydraulic dental chairs are world renowned for providing reliable long-term performance for operators, as well the smooth comfortable seat motion that minimizes patient stress and anxiety as they are gently glided into treatment position.



Dental Light -Realization of Natural Light

LED dental lights produce natural light. By placing LED bulbs in round pattern, our LED lights produces natural light with high color rendering and performance, for correct tooth shade matching and for differentiation of soft tissue.



Dental X-ray -Safety and Reliablity

Our reliable quality is based on the strict inspection standard and user-friendly designing such as our easy to understand operation panel, smooth arm movement and easy to position compact tube head. These unique features are helping us gain a worldwide reputation.

Since 1921 in Osaka, Japan, Takara Belmont has been casting metals for chairs.

In the future, we will continue our ongoing efforts to actively introduce advanced manufacturing facilities and efficient production systems in response to the changes in times, produce high-quality, innovative products and enhance the reliability and safety of our products.









Upholstery Processing Technique

In responding to the customer's demand, highly trained upholstery workers handcraft the chair by using staple gun and pressing machine.

Production Inspection

To provide high quality products to all customers, we provide inspection of electricity and functions to the products before shipping.

Machining Technology

We process a broad range of components from cylinders, ram shafts, to many precise components by ourselves. Core components are inspected for dimensional accuracy at micro-meter level.



Committing to contribute to the healthy life not only for patients, but also for dentists & their team.

Ergonomics in Dental Practice

Dentists are at high risk for musculoskeletal disorder due to the nature of their work. 81 percent of dentists surveyed said they had experienced work-related neck pain during the year, 51 percent reported lower-back problems, 44 percent said they had shoulder pain, 43 percent had upper-back pain, 38 percent reported hand pain, 30 percent mentioned mid-back pain and 14 percent indicated arm pain.

The key to preventing work-related musculoskeletal disorder is ergonomics - the science of fitting the work environment to the worker.





Patient-friendly Design

What the dentistry can offer is to propose a lifestyle, in which people can enjoy a healthy oral environment throughout their life.

We, Takara Belmont, are always ready to support dentists and all the staff by designing the space and coordinating dental equipment to invite not only patients but also people, who have healthy teeth to learn about preventive oral care.

For example, our leg-folding chair. The patients can easily take a seat as if they sit on chairs at home. It helps to create a more comfortable space for the patients. Takara Belmont is the first manufacturer of the leg-folding chair in Japan.

As a company advocating the value of "Health and Beauty", we respond to people's demand by designing space, developing products and services, and hope to support them to maintain the happy teeth.







eurus S8



EURUS

With a name inspired by the East Wind, which brought blessed rain to grain fields in ancient Greek mythology, Eurus embodies our belief of bringing highly innovative, reliable products to the world of dentistry.

Compact Space

Wide Variety of Treatment Styles

The ambidextrous dental unit is ideal for both right and left handed dentists and assistants. The changeover is easy and takes only a few minutes. This function makes the unit suitable for a variety of operating styles and increases efficiency. For example, the rotatable assistant holder tray offers dentists a more comfortable operating space and allows for two-handed operations. Positioning the doctor table behind the chair enables assistants to prepare during consultations.





Two-handed position
The long reach and light movement of
the arm facilitate two-handed treatment
to be administered without stress.





Light As Air

A highly manoeuvrable and lightweight doctor table reduces fatigue for dentists. It allows the doctor table to be moved with only a minimum of effort. This characteristic can lead to a prolonging of a dentist's career. The doctor table and cuspidor can be moved over a wide range from left to right. Belmont has over 25 years experience in developing this wide range of motion specification.



Non-standard holder (indent)

Smooth, Comfortable and Relaxing

Chair

The advanced hydraulic chair mechanism ensures the smooth, quiet and efficient movement of every chair function. This allows the patient to be manoeuvred effortlessly into the most comfortable position to receive their treatment.

Accommodating the unique age, height and stature of every patient, from small children to adults, the position of the new twin-axis headrest system can be fine-tuned in increments of 7.5 degrees for the ultimate in precision and comfort.

Featuring new advancements in manoeuvrability, the headrest can also be height adjusted effortlessly by the dentist and assistant to ensure complete comfort for the patient.

The chair has a low initial height with 200kg load capacity.



Rotating armrests

The armrests are designed with a length and shape that make it easy for patients to hold, and the side-rotating mechanism helps elderly patients mount and dismount the chair with ease.



Wireless foot controller with chair stick switch Option
The dental chair can be entirely manipulated in a hands-free, hygienic and intuitive manner using only the operator's foot.



Power headrest Option

The full 120mm stroke and the neck angle can be easily adjusted without imposing excessive strain on the patient. The motorised headrest can also be manoeuvred hands-free via the foot switches, which ensures complete hygiene and safety.

Fit The Height

Arm Holder

The table's vertical range of motion is 840, 865, and 890 mm, adjustable in three-stage height.



Standard holder

Handle - ideally suited for touch-panel operation The handle is not only designed for an easy hold-and-move motion, it also provides an ideal resting position for the hand. This reduces stress and fatigue, and it facilitates reliable touch-panel operation.





Instrument holder's position

Designed to be independent from the table, the instrument holder can be freely adjusted via a three-axis adjustment mechanism to achieve the best working position for the dentist



Non-standard holder (indent)

Hose guide

The smooth, frictionless motion of the roller frees the dentist and dental assistant from the weight of the hose when drawing out dental instruments. This simple, yet ingenious design minimises strain from a working process that is repeated frequently throughout the day.

Your style at the touch of a finger

Interface

Designed to improve efficiency and reduce stress during treatments, the screen only displays necessary information when a handpiece is selected and picked up. This allows the dentist to concentrate on the treatment at hand as well as the patient, without unnecessary distraction.



Motor

Configuring the screen for handpiace motors, which is capable of storing four different settings, allows for Intuitive, seamless operation. RPM, torque and rotational direction can also be selected with a single touch of the screen.



Sideways operation

A new slide foot controller allows for precision adjustment via sideways position motions and step-on activation. Maximum RPM can be set for the handpiece quickly and intuitively.

Providing Quality

Cuspidor

The cuspidor has been expanded to offer larger capacity and easier access which, with the separator, ensures the highest levels of hygiene are maintained.









LED Indicator

The unit features a safety Indicator that displays the lock-status of the dental chair and handpieces in an easy-to-see way. This ensures the dentist and dental assistant quickly and easily understand the status of the unit.

Suction system

The unit provides effective and reliable suction that has been significantly increased by 185%. This system can also be adapted for twin suction, which proves useful for more intensive surgical procedures.

*Comparison of our conventional products under our conditions.

Improved cleaning performance

Cleanlines

Hygiene inside the chair unit is maintained by water line flushing, which means that patient-friendly treatments can be administered in complete safety and confidence. As an environmentally friendly option, the unit comes with a separator, or a separator with an amalgam collector, to reduce the need for everyday maintenance.





Autoclavable holder
The dentist outlets and the assistant outlets can both be autoclaved and detached for cleaning.

Suction line clean system Non-standard (indent)
Daily use of the integrated suction
cleaning system prevents build-up of
debris in the suction line and helps to
maintain constant suction.





action for easy cleaning.



Detachable bowl The ceramic bowl is easily detached for cleaning to make sure it is always clean and hygienic.



Water bottle Option
Clean water is easily accessible for handpieces and the cupfiller via a water bottle that's conveniently located and secured immediately under the doctor table.

Bright and gentle

Lighting Fixture

Our lighting fixtures are designed to reproduce natural light with a high functional colour rendering capacity. The light's composite mode can be active, and the light can be turned on and off, via its touch-less sensor.



LIGHT UP BRILLIANTLY

The new 900 dental light provides variable intensity up to 39,000 luxes with rectangular light pattern and is bright enough to operate in the latest treatment styles.

It is easy to be turned on and off by the touchless sensor, and the composite safe mode can be active as well.

Flexible lighting

The Belmont spring arm allows the light to be moved around effortlessly and remain in the desired position, so that the dentist or dental assistant can begin treatments without delay. The arm spans a wide area, including the maxillary and mandibular areas of the patient, to ensure optimum illumination for every dental procedure.



Colour variations

Colour coordination

The upholstery can be colour coordinated to suit the interior aspects and design of your dental clinic from four themes: Warm, Cool, Earth and Vitamin; and from two upholstery styles: seamless and luxury.

[Seamless]

Warm Colour

Coral Pink /ES20 Light Pink /ES30

Beige/ ES21

Camel /ES14 Bordeaux /ES 12 Black/ ES8



Cool Colour

Mint Blue /ES5

Indigo Blue /ES25 Light Blue /ES16

Purple /ES26 Light Gray /ES33 Gray Teal /ES28



Earth Colour

Rose /ES34 Mauve /ES27 Forest /ES31

Aldo Blue /ES23

Ash Gray /ES32 Chocolate /ES 11



Vitamin Colour

Scarlet Red /ES24 Mandarin Orange /ES22 Lemon Yellow /ES2

Yellow Green /FS18 Sky Blue /ES7

Blue /ES29



[Luxury]

Black /SH1 Dark Blue /SH9 Maroon /SH7 England Green /SH18

Orange /SH22 Camel /SH10 Mustard Yellow /SH6

lvory /SH15

Select a Eurus unit that's best suited to the specific needs of your clinic and patients from a wide range of options.



Variations

Conceptualised by Belmont to capture essences of cleanliness, tenderness and supreme ergonomics, the Eurus Series presents a comprehensive range of stools, dental chair units and dentist tables to accommodate the individual requirements of you and your dental team as well as patients.

With a proud and distinguished heritage in the dental industry, Belmont prides itself on creating innovative products that are conceived through a deep understanding of the professional and treatment requirements of dentists and dental assistants, and the expectations of your patients. Our appreciation of established and emerging technical demands and dental procedures, and the anxiety experienced by many patients, are also reflected and inherent in every equipment concept to make the delivery of treatments, and a visit to the dentist, as pleasurable and comfortable as possible.

S8 Standard

[Chair] Contour Chair

[Unit Mount Type] Rotation Cuspidor & Base Mount Dr. Table Rotation Cuspidor & Cart

[Delivery Type] Holder/Place

[Light] Unit Mount LED Light



S8 VacPac

[Chair] Contour Chair

[Unit Mount Type] VacPac & Base Mount Dr. Table VacPac & Cart

[Delivery Type] Holder/Place

[Light] Unit Mount LED Light





S8 Without Cuspidor

[Chair] Contour Chair

[Unit Mount Type] Base Mount Dr. Table Cart

[Delivery Type] Holder/Place

[Light] Unit Mount LED Light



clesta elli swing



Basic, Light, and Flexible

Basic

- Easy changeover to right/left handed use
- Easy-access handpiece holders
- Compact space
- Smooth tilting movement by hydraulic system
- Ergonomically shaped cushions

Light



- Reducing fatigue for arms of dentists and assistants and leads to prolonged working lifetime

- Fitting position with precise movement by light headrest

Flexible

- Making comfortable operating space
- Providing a variety of operating styles and preparation zones as you like



Light as Air

The doctor table is simple and practical to use. It can be moved lightly and easily.



Just "SWING it" to Right/Left Handed Use

This model provides two- and four-handed working positions and various operation and preparation zones.

Ambidextrous changeover is easy and takes only a few minutes.

Dentists and assistants can position the doctor table and assistant holder as they like.





Fit your position
The doctor table can be adjusted with three-level height and each dentist can position it appropriately.

Smooth and Relax

The hydraulic mechanism ensures smooth and quiet movement of the joints and provides a comfortable feeling for patients. The new twin-axis headrest system, capable of making a wide range of adjustments, can be fine-tuned in increments of 7.5 degrees, accommodating patients from children to adults. The chair has a low initial height with 200kg load capacity.





Rotating armrests

The armrests are designed with a length and shape that make it easy for patients to hold, and the side-rotating mechanism helps elderly patients mount and dismount the chair with ease.

Dental Light



EURUS LIGHT

The EURUS LIGHT provides sensor activation, stepless intensity between 3,000 lx and 28,000 lx, reproduces natural light with a high functional color rendering capacity.





Detachable light components

Lens cover and autoclavable light handle can be removed with a single lever action for easy cleaning.



Specifications

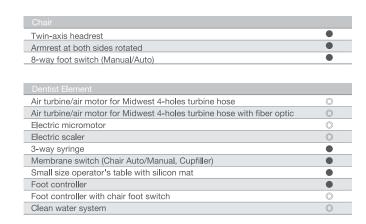

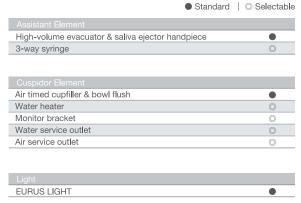



Upholstery Colour Options


The colors come in four themes: warm, cool, earth and vitamin. Any choices from these leather colors should match the interior appointments of your clinic.



Dimensions

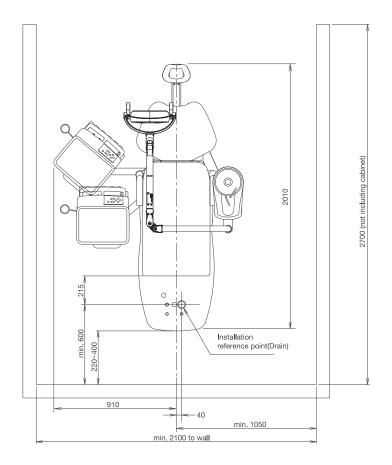

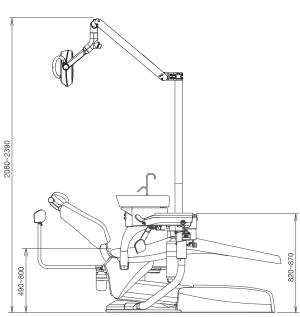

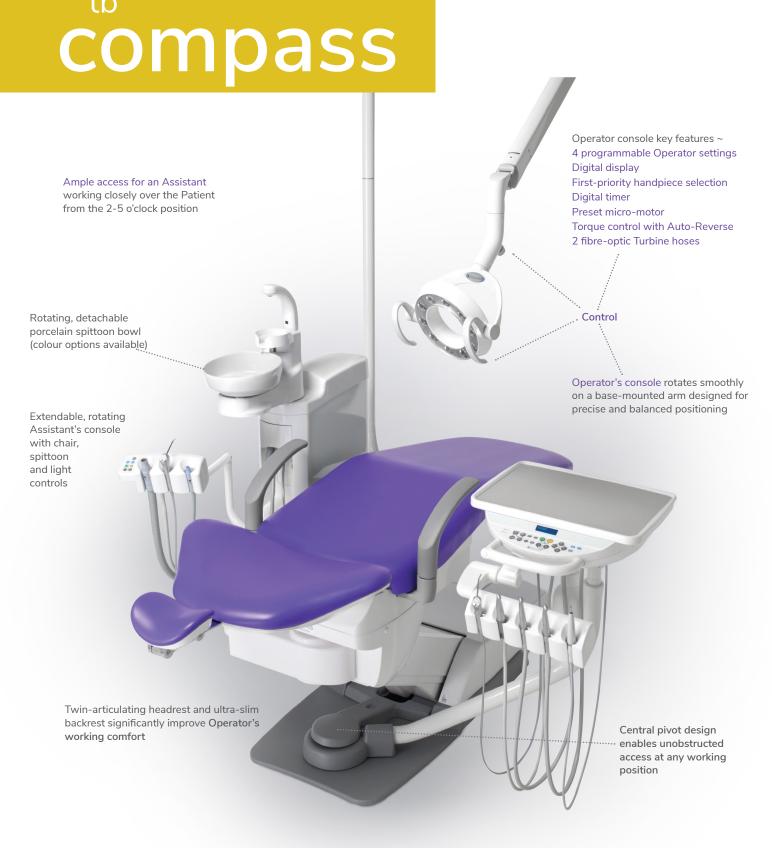

Control and discretion are the finest traits of the tbCompass delivery system. An intelligent delivery solution which seeks to conceal the instruments from Patientview without comprising working comfort. Designed for uninterrupted dentistry in the most demanding surgery environments, tbCompass is pointing to the future.

Ambidextrous delivery options for better daily dentistry

Unrivalled console mobility and tbCompass's determined versatility bring all the benefits of convenience, flexibility, efficiency and Patient comfort.

Belmont's base-mounted delivery system design is an understated advancement in dentistry and Patient care.

Ambidextrous capability across the range doubles the qualification of tbCompass in any configuration to perform excellent service for any Operator and procedure.

Key Features

- Ambidextrous
- Below-the-Patient swing-arm delivery system
- Ideal console arrangement for Assistant's 'preparation zone'
- Vac-Pac Systems option
- Mobile delivery systems option
- Air and Electric operation system
- Seamless or Plush upholstery option
- Standard or Ultrasoft upholstery colour options

tbCompass 'Surgery System' (Ambidextrous)



tbCompass 'Flexible Dr System' (Ambidextrous)

chair and Operator's console



tbCompass 'Flexible Cuspidor System' (Ambidextrous)

chair and Operator's console



compass

Complete package (unit, patient chair, assistants holder & LED overhead light)

Package Inclusions

- Central HVE & SE with cut off valve
- 2 x DCl 3-way syringe provision
- 2 x fibre optic tubing
- 1 x NSK NLX nano provision
- 1 x NSK varios 170 scaler provision
- Installation 2.0L bottle
- ratchet gear arm
- service outlet for water & air
- male connector for water & air
- swing arm (w/monitor hole, w/out cap)

tbCOMPASS chair

- Twin axis headrest
- Left & right armrest with foot switch
- Colour: 2/ES8 black, 2/ES8 black



compass

unit, chair and assistants' holder

Package Inclusions

- Dental unit & chair tbCOMPASS 230V/50HZ for flexible DR system without cuspidor unit
- Holder type handpiece
- 1 x DCl 3-way syringe provision
- 2 x fibre optic tubing
- 1 x NSK NLX nano provision
- 1 x NSK varios 170 scaler provision
- Installation 2.0L bottle

tbCOMPASS chair

- Twin axis headrest
- Left & right armrest with foot switch
- Colour: TBA, Other specifications are as per factory default
- D.C.I 3way syringe (3439)
- Cover hose
- Spacer



clesta elli



variations




Holder type Rod type



Cart Holder type



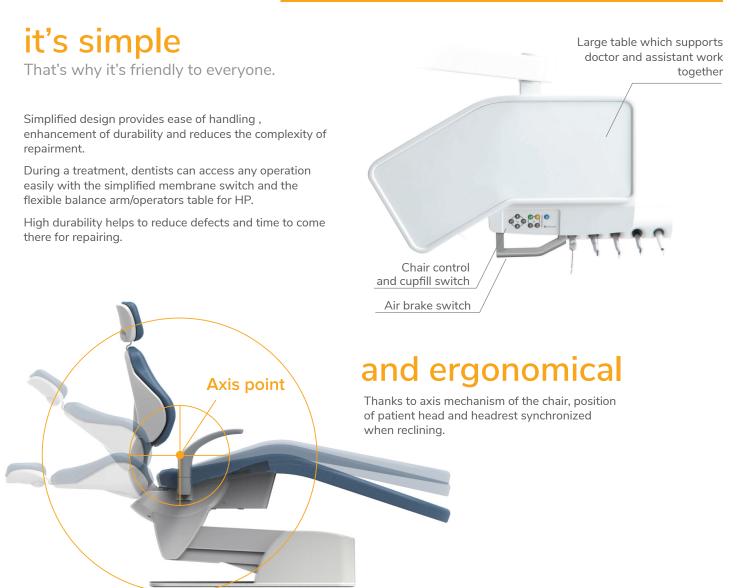

34

options and accessories



Flushing tray

The flushing system enables easy cleaning of handpiece water lines. The flushing tray as an extra option allows hands-free flushing.



Monitor bracket

The flat screen monitor bracket can be adjusted to suit individual need. (not including monitor)



Clean water system

The clean water system offers a choice of 750-ml or 2-litre bottles and can be integrated into the unit to supply handpiece water. The aesthetic cover is an extra option.



Rotating armrest

A swing-out armrest for the right side makes it easy for patients to hold and helps elderly patients mount and dismount the chair with ease.



Assistant's console

Assistant sub-tray increases assistant's work efficiency. (Metal sub-tray is not included.) Moreover, curing light holder enhances the operability of the curing light. HVE tips and SE tip are optional.



Double articulating headrest

Suitable for children and adults, and it's improved to adjust patient positioning easily.

unit and chair

Cuspidor

Easy-to-clean porcelain cuspidor bowl is shaped to prevent rinsing water from spilling out.



- Powerful hydraulic system enables a wide range of lift heights.

Balance arm

Employing an air brake system, the balance arm enhances mobility and flexibility during dental practice.

Operator's table

Large operator's table ensure ease of use. Swivelling hardpiece holder for most appropriate position for dentist.

- Ergonomically designed chair cushion ensures patient comfort. The reclining backrest is designed to prevent the patient from sliding.

light



300 LED LIGHT

The 300 LED light provides sensor activation, stepless intensity between 3100 lx and 28000 lx, three axis rotation and clear light pattern



Flexible Lighting

The Belmont spring arm allows the light to be moved around and stay exactly where it is desired so that the dentist/dental assistant can start treatment smoothly. The arm covers a wide area, including the maxiliary and madibular areas of the patient, for sufficient illumination.





Monitor bracket

The flat screen monitor bracket can be adjusted to suit individual need. (not including monitor)



clesta II



variations

Chair Mount -original-

Over-patient versatility for comfortable singleand four-handed working positions. Ergonomic and hygienic features are standard on all CLESTA II treatment centres.

Shown as

-Unit : CLESTA II

-Chair : CLESTA II CHAIR -Light : 900 LED LIGHT



Axis backrest

Synchronized axis backrest ensures optimum patient comfort at all times. Ultra slim backrest for ease of access for both operator and assistant.



Durable diecast base

The sleek styling of the durable aluminium diecast base offers the dentist improved access with maximum chair stability.

Pedestal -CLAIR-

Pedestal type cuspidor unit has no umbilical; this eliminates the need for a utility box, contributing to a tidier appearance and easier cleaning.

Shown as

-Unit : CLESTA II -Chair : CLAIR CHAIR









Pedestal type / Cart



Adjustable chair positioning

A wider range of movements offers a practical working environment that meets all the operators' requirements.



Service centre

The built in service centre eliminates the use of external umbilical hoses.

Intuitive Control - technology you can rely on

E Type Electric Type Control Panel



- 1 Power Indicator
- 2 Safety Lock Indicator
- 3 Coolant Water ON/OFF Switch
- 4 Light Pack Switch
- 5 Micromotor Speed Setting Switch
- 6 Micromotor Direction Control Switch
- 7 Function Switch
- 8 Store Switch

- 9 Decreasing Switch
- 10 Increasing Switch
- 11 Chair Auto Mode Switches
- 12 Chair Manual Control Switches
- 13 Cupfiller Switch
- 14 Bowl Flush Switch
- 15 Dental Light Switch
- 16 Function Indicator

Function Switch (E Type)



Timer

Timer set for alarm.



Group selection mode

Preset position for chair and preset rotation speed for micromotor can be set for 4 individual dentist.



Flushing system

The E type is available with a short and long time flushing system.



Control panel switching sound ON / OFF

Sound can be deactivated.



Fiber optic handpiece lighting mode

Handpiece lighting mode can be chosen.

Digital Indicator (E Type)




Micromotor mode

The micromotor rotation speed can be selected by sliding the foot controller pedal right or left. Rotation speed is displayed.



Electronic sound for timer

4 electronic timer sounds to choose from.

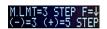

Micromotor maximum speed setting (optional)

The maximum rotation speed of the micromotor can be selected in 3 steps or 5 steps.



Coolant water ON / OFF switch

When using micromotor(s), either 2-mode or 4-mode of coolant water can be selected by mode select setup.



Cupfiller and bowl flush

Select cupfiller and bowl flush to operate together or independently.



Scaler mode

Easy to confirm which mode of scaling has been selected (scaling, end, perio). Power output is displayed.

Upholstery Colour Options

The colors come in four themes: warm, cool, earth and vitamin. Any choices from these leather colors should match the interior appointments of your clinic.



Light



900 LED Light

The 900 LED light provides sensor activation, steppless intensity between 4000 lx and 35500 lx, fluid three axis rotation and clear light pattern with high color rendering.

Options & Accessories



Flush out tray



Flat screen monitor bracket (not including monitor)

Instrument Delivery







Unit



Integrated cuspidor

The highly durable integrated cuspidor creates an attractive, aseptic environment for the patient and the smooth surface makes it easy to clean. Detachable glass bowl is durable, transparent and easy to clean.



Operator's table

Incorporating holder type instrumentation with "first priority system"*, the operator's table is adjustable both vertically and horizontally to give optimal positioning.



Assistant's console

The assistant's console includes 3-way syringe, vacuum and saliva ejector tubings, and easy clean membrane switches to control chair movements, cuspidor bowl flush and cupfiller.





Foot control

A2(control of air instruments, water on/off and chip air) for A type. SE(control of air instruments, water on/off, chip air and control of electric motor/scaler) for E type. Foot controller with handle is available as an option.


Dimensions


Chair Mount -original-

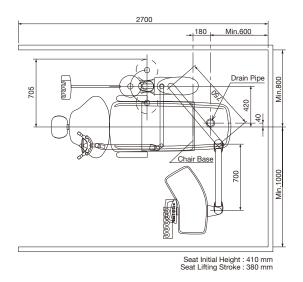

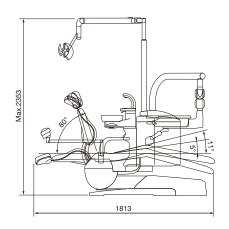


Pedestal -CLAIR-

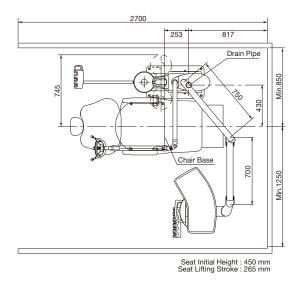


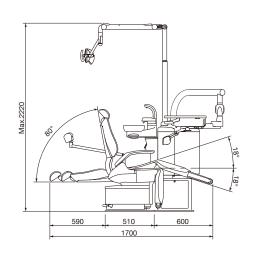

Rod Type

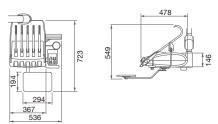

Cart Type

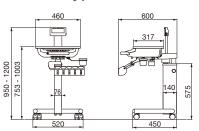

Names of Components

Dental Unit : CLESTA II

: CELEB

Dental Chair : CLESTA II CHAIR
: CELEB CHAIR

: CLESTA II CHAIR (EURUS TYPE) : CELEB CHAIR (EURUS TYPE)

: CLAIR CHAIR

Dental Light : EURUS LIGHT

: 900 DENTAL LIGHT TYPE 920 : 900 DENTAL LIGHT TYPE AL-920

: 300LED DENTAL LIGHT TYPE 320S : DENTAL LIGHT MODEL 720S

eurus

swivel chair

EURUS SWIVEL CHAIR allows patients to enter and exit the chair with an experience similar to sitting in a home armchair. The unique design is ideal for face-to-face consultation and communication between dentist and patient, permits direct frontal access to patients for treatment procedures.



With the legrest folded away, dentist and patients can interact more comfortably during face-to-face consultation and treatment procedures that may require direct frontal access, such as in cosmetic and orthodontic dentistry and for the assessment of jaw alignment.



The chair rotates on its base left-to-right so patients can be repositioned for taking X-rays for treatment procedures. The chair's small footprint in the upright position requires less space within the surgery.







Synchronised Seat Motion

Synchronised shockless movement of the backrest, seat and legrest enhance patient comfort as they are positioned for treatment. As the backrest reclines, the seat tilts backward and the legrest slides outward to provide support for the patient's legs.

Chair Rotation

The chair rotates 105 degrees to the left & right on each side of center, and patients can be repositioned as needed so that treatment and consultation can be done at an appropriate position.





Armrest

Armrest rotates outward to provide excellent access for patients when entering or exiting the chair. With the chair base bolted to the floor for maximum stability, armrest can provide support for patients when being seated and when getting up from chair.

Headrest Choices


The double articulating headrest version is manually adjustable and provides a wide range of positioning flexibility for superb access to patient maxillary and mandibular arches for examination and treatment. The power headrest version permits hands-free positioning of the patient's head for examination and treatment.




Seat Positioning Controls and Rotation Lock

Foot activated stick switches are integrated into the doctor and assistant side of the chair base for convenience and ease of seat positioning, without the need to resort to a wired external foot control. A pivoting foot pedal located on the side of the chair base is utilised for hands-free locking/ releasing of chair rotation.

dimensions

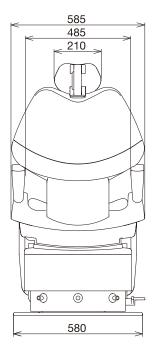

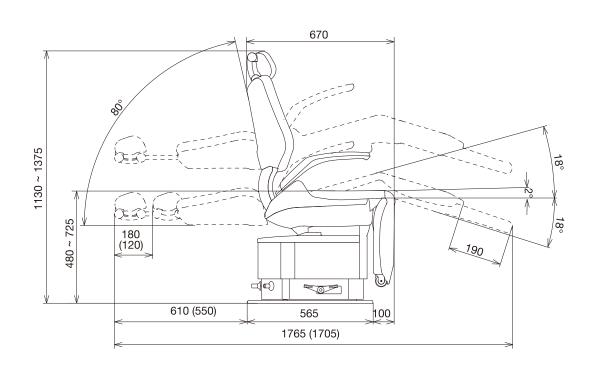

46

Knee-break -Clair-



Folding legrest chair

Slim, yet comfortable backrest and chair enable easy access to the patients.



Armrest

Swing-out armrest allows the elderly, handicapped, and wheelchair patients to enter and exit the chair safely and with ease.



Chair stick switch

The stick switches enable the operator to position the patient quickly and accurately, keeping hands-free and clean.



The ergonomic, axis mechanism backrest supports the patient comfortably and accommodates the natural movements of the patient's body.




a spectrum of

upholstery colours

to set your dentistry mood



Warm Colours

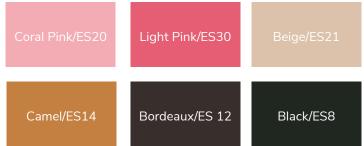

Cool Colours







Earth Colours

Rose/ES34

Mauve/ES27

Forest/ES31

Aldo Blue/ES23

Ash Gray/ES32

Chocolate/ES 11

Vitamin Colours

Scarlet Red/ES24

Mandarın
Orange/ES22

Lemon Yellow/ES2

Yellow Green/ES18

Sky Blue/ES7

Blue/ES29



dental light 900



900 DENTAL LIGHT

provides cool white LED illumination of the oral cavity for correct tooth shade matching and for the differentiation of soft tissue structure.

The well-defined, 85mm x 155mm rectangular light pattern remains unaffected when an obstruction of up to 50% of light head is placed between the light and patient.

900 DENTAL LIGHT operation is sensor activated so it is not necessary to look away from patient to operate light. Touchless ON/OFF and Composite Safe Mode activation are done without ever contacting the light.





Long Life

With LED bulb life expectancy of 40,000 hours, it is likely that bulb replacement will never be needed.

Economical

Low power consumption and no heat generation add to long term savings and to patient comfort.

Cleaning

Clear bulb cover ring prevents mists and dust from getting inside of the light head. The entire light head is easily cleaned.

features

Brightness distribution

900 DENTAL LIGHT is engineered to minimize eye fatigue. The rectangular light pattern is brightest inside of the oral cavity and gradually becomes less bright across the patient's cheek area.



Fluid positioning

Smooth, three-axis light head rotation permits stress-free light head positioning with one hand. Positioning flexibility directs illumination to area where it is needed.



Intensity adjustment

Each operatory has different ambient light conditions. With this in mind, 900 DENTAL LIGHT lights include a stepless intensity control to adjust brightness to fit the actual operatory light condition. Light output can be adjusted between 4,000 lx and 35,500 lx, while maintaining the 5,000°K colour temperature that is ideal for shade matching.



installation variations








(post not included)

dimensions

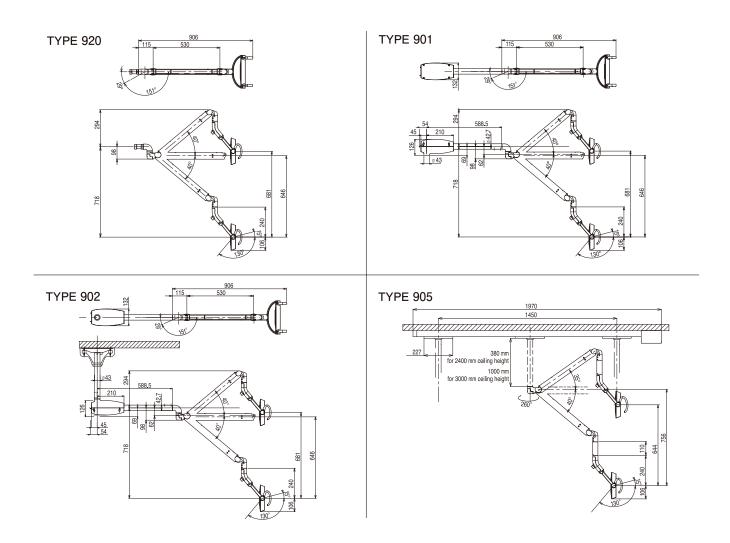

specifications

Rating Voltage	100~240V 50/60Hz (AL920:DC20V)	Colour rendering index	More than 90
Intensity	4,000 ~ 35,500 lx (Stepless)	Focal Distance	650 mm
	5,500 lx at composite safe mode	Light pattern	85mm x155mm
Colour temperature	5,000° Kelvin at any intensity	Bulb Type	High Rendering LED

The manufacturers reserve the right to alter and improve the specification without notice.

Phot-XIIS



Premium

LCD touchscreen sub-controller



Intuitive and easy selection

Intuitive and easy selection of tooth type, cone type, film / digital modes, film speed, kV, mA and patients size.

Customizable LCD display view

LCD display view can be customized for operator preference. Built-in USB side port enables operators to import personalized screen saver photographic images or dental practice / doctor's name for practice identification and branding.

(A) Sum of radiation dose display


Check the radiation dose of each patient. Reset the sum by keeping touch the sum number.

B Display of exposure time


Change the exposure time easily by touching the arrow and obtain the image you desire.

C Display of radiation dose of each exposure.


Sensitivity setting

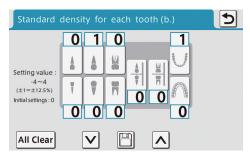

Select your image receptor on the screen, then suitable sensitivity setting will be selected from 16 steps automatically.

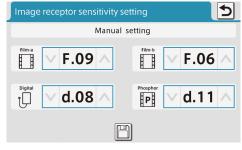

If the darker (or lighter) radiographs are preferred, the sensitivity setting can be adjusted manually

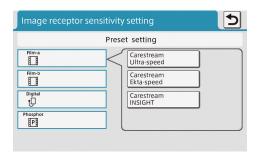

Tooth density

Tooth density ratio settings can be adjusted from the default setting for each tooth

Reliable DC X-ray by Belmont

Belmont's High Frequency Technology allows for


- · High resolution image by a small focal spot
- Reduced skin dosage with optimum radiation level for every diagnostic need
- Easy to position compact tube head
- · Ideal matching for digital imaging system






Features

Small 0.4 focal spot value for the sharpest image possible

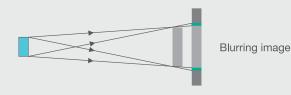

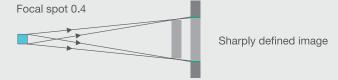

Radiation dosage

 $\mbox{\sc PHOT-XII}\mbox{\sc S}$ utilizes a Direct Current (DC) generator to minimize radiation absorbed by patients.

In addition, there are 36 programmable exposure time settings between 0.01 to 2.0 seconds to fine tune the precise amount of radiation required to produce an excellent radiographic image.

Easy to positon compact tube



The tube head is compact and ear This enables the the tube head as

Realization of slim and light a



Our spring arm is designed to have a good balance any position within the movable range and make if at any point when the hands are released. Moreoveleps to reduce the vibration of the tube head after These methods realize the smooth and easy positions.

Standard Sub-controller





head

s lightweight, sy to hold. operator to align ccurately every time.

m



Light head, smooth movement

ce with load and repulsive force at t possible to make tube head stay ver, well balanced rigidity design er positioning. ioning of the tube head.

Optimum radiation level for every patient

The operation is simple. Select patient type and tooth type, then press the exposure switch. The exposure time can also be adjusted manually.

After the exposure, the radiation output can be displayed in the window.

Tube voltage selection

PHOT-XIIs provides 2 tube voltage settings to optimize imaging for specific purposes:

[60 kV Setting]

- Setting provides balanced contrast for general caries diagnostic imaging.
- Setting provides excellent definition and detail for periapical and bone structure diagnostic imaging.

[70 kV Setting]

- Setting provides excellent diagnostic imaging for periodontal disease examination.
- Setting provides a wider latitude and longer gray scale.

Five types of installations

- WK Wall mounted type
- RK Room type with stool
- FM Mobile type with castors
- FK Floor mounted. Sub-controller can be installed on the pole (FK1) or on the wall (FK2)
- CK Ceiling mounted type. Dental light can be mounted on the pole as an option (CKL Type)





- WK Type (Sub-controller can be mounted on the main controller or on the wall)

RK Type







FM Type

FK1 Type

CK Type

Options


- Horizontal arm length of WK Type is selectable: 300mm, 500mm, 650mm, 800mm (standard) and 1000mm
- Ceiling pole length of CK / CKL Type is selectable: $600 \sim 1350 \text{mm}$
- Long cone (Short cone + Long cone collimator) f or paralleling technique
- Rectangular collimator
- Hand exposure switch with coiled cord (Option for WK, RK, FK2 and CK. Standard for FM and FK1.)











dimensions

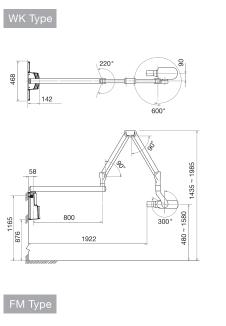

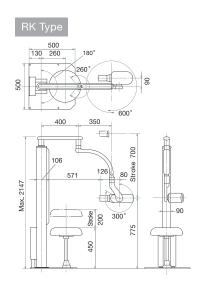

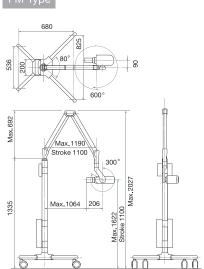

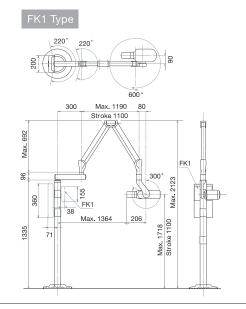

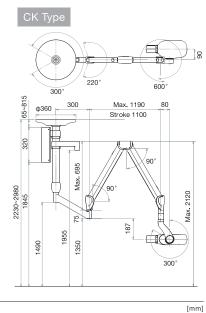

ſ	r	۲	١	r

Nominal focal spot value	0.4 (IEC 60336)	
Rated tube potential	60 kV / 70 kV (selectable)	
Rated tube current	3 mA / 6 mA (selectable)	
Rated line voltage	100 / 110 / 120 / 220 / 230 / 240 V AC 50 / 60 Hz	
Rated line current (at 70 kV / 6 mA)	·11.0 A (100 V) ·10.5 A (110 V) ·10.0 A (120 V) ·6.4 A (220 V) ·6.2 A (230 V) ·6.0 A (240 V)	
Line voltage range	±10%	
Exposure time	0.01 to 2.0 sec.	
Total filtration	2.0 mmAl Equivalent at 70 kV	
	Regular cone	203 mm
Source to skin distance (SSD)	Long cone (option)	305 mm
	Rectangular collimator (option)	SSD of cone + 40 mm
Duty cycle	1:30	
Half value layer	1.5 mmAl over	
Donoity proprogram	Film	16 programs
Density preprogram	Digital Sensor	16 programs

Names of Components Dental X-ray :PHOT-XIIS MODEL 505

:PHOT-XIIS 505 :PHOT-XIIs MODEL 505 :PHOT-XIIs 505 :PHOT-XIIs 505LCD

:PHOT-XIIs MODEL 505 WITH ACCESSORIES

Stools



Doctor & Assistant Stool Package

Comfort

Meets the needs of various treatment styles and supports healthy posture.



Doctor stool EURUS STOOL

- Seat height and tilt adjustment by gas cylinder.
- Bucket-shaped seat and backrest help treatment process over a long period.
- Tilt functions for backrest and seat are available. (Backrest and seat tilt separately.)
- Three casters with brakes prevent unexpected move.

[SPECIFICATIONS]

Model	AH-012D
Weight	12 kg
Initial height/stroke	440/155 mm
Tilt control of the seat and backrest	Manual operation



Assistant stool DH-007N

- Seat height adjustment by gas cylinder.
- Smooth and light touch operation for seat height adjustment by hand lever.
- Simple and easy mechanism makes added flexibility possible.

[SPECIFICATIONS]

Model	DH-007N
Weight	10 kg
Initial height/stroke	465/195 mm
Armrest/backrest rotation	360° by Manual operation

Doctor & Assistant Stool Package (Cont)



Backrest height choice

The height of backrest can be selected from 3 positions, 30 mm difference between each position, to support your healthy posture.







Casters with brakes

Brakes are available for 3 of the 5 casters to prevent the stool from moving unexpectedly when no load applied, which allows the doctor sit on the stool safely. When you sit on, the brakes will be unlocked. When you stand up, the brakes will be locked.

Tilt control

The seat and backrest tilt separately to maintain the ergonomic sitting posture. Seat tilt can be adjusted from horizontal to 9-degrees forward, sitting in ideal position to minimize pressure on thigh and support posture.

Color variations

The colors come in four themes: warm, cool, earth and vitamin. Any choices from these leather colors should match the interior appointments of your clinic.




[Luxury]





Dimensions

DH-007N

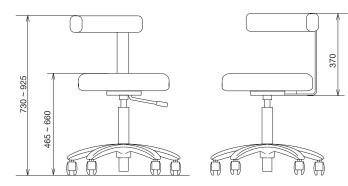

LS-007GSB

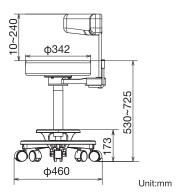

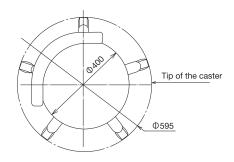

Doctor Surgical Stool

Comfort

Meets the needs of various treatment styles and supports healthy posture.



Doctor stool LS-007GSB

[SPECIFICATIONS]

Model	LS-007GSB
Weight	13 kg
Height	530mm-725mm
Base	
Backrest	Standard



Rise into beauty. Walk into life 100 years of Takara Belmont

Across all divisions of TAKARA BELMONT, we have tried to accurately read and reflect local market needs in new equipment.

Through our worldwide network of offices and subsidiaries, Takara Belmont continually monitors current industry and market needs as we strive to design and introduce high quality, creative products that help increase work efficiency and productivity.

DISTRIBUTORS

Contact One Of Our Trusted Dealers Below For A Great Deal

NSW

Mark's Dental Pty Ltd Mark Mangesi P 0404 448 616

Urdent Dental Equipment Services Jamie Warner P 0411 217 278

Australian Imaging Sean Thomson P 0428 909 539 Daniel Hatfield P 0438 747 666 www.australianimaging.com.au

NSW - NORTHERN & NORTH COAST

Country Dental & Autoclave Services Ben Hargraves P 0418 657 203 www.countrydental.com.au

NSW - HUNTER REGION

Hunter Dental Tech Services Owen Kurtz P 02 4946 8007 www.hdts.com.au

VIC

Borg Dental John Borg P 03 9329 6276 www.borgdental.com.au

Medi-Dent P 1300 886 674 medident@medident.com.au

WA

Med and Dent Luke Ferris P 0411 592 821 David Jones P 0475 693 242 www.medanddent.com.au admin@medanddent.com.au

QLD

Leading Dental Steve Symonds P 07 3139 1549 www.leadingdental.com.au

TAS - (ALL AREAS)

Medi-Dent P: 1300 886 674 medident@medident.com.au

NEW ZEALAND Complete Healthcare Solutions Praneesh Prasad P +64 2171 5107 www.completehsl.com

